

ETERNITY

POWERFUL RUBY LASER FOR RESISTANT INKS



ASTANZA

SERVICE | TRAINING | MARKETING SUPPORT | TECHNOLOGY

ETERNITY^{TSR}

POWERFUL RUBY LASER FOR RESISTANT INKS



Engineered for power and reliability, the Eternity^{TSR} takes tattoo removal practices to the next level.

DESIGNED FOR TARGETING STUBBORN INKS

About 20-30% of tattoos feature colors such as lime green or sky blue or have faded from black to green over time and pose a difficulty for laser tattoo removal. The Astanza Eternity^{TSR} precisely targets these colors – in addition to other resistant shades of blue, green, violet, and black. See fast, effective removal on the most stubborn inks with the Astanza Eternity^{TSR}.

POWERFUL RUBY WAVELENGTH

A laser's wavelength is one of the most important factors for targeting tattoo ink. The 694 nm wavelength of the Eternity^{TSR}'s ruby laser is excellent at treating green, blue, violet, and black inks that other lasers leave behind.



HIGH ENERGY OUTPUT

Boasting a full 1.2 J of energy, a 2 Hz repetition rate, and up to 20 J/cm² of fluence, the Eternity^{TSR} is the most powerful ruby laser on the market – quickly removing the most stubborn ink colors with ease. The Eternity^{TSR} has up to 6x the power of dye handpieces and requires no expensive consumables.

IDEAL FOR RESISTANT COLORS



Before



After 9 Tx with competitor system
at 1064/532/650/585 nm



After 3 Tx with Eternity^{TSR}*

**Images courtesy of Erase the Ink M.D.*


IMPRESSIVE RESULTS FOR ELITE PRACTICES




Before



After 6 Tx**

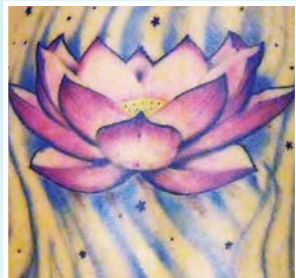

Before

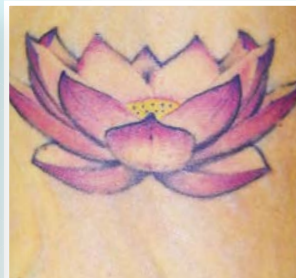

After 6 Tx**



Before



After 4 Tx**



Before



After 1 Tx**

The Eternity^{TSR} is also highly effective for the treatment of pigmented lesions.

CUTTING-EDGE FEATURES

The Eternity^{TSR} combines a variety of advanced features as the world's most powerful ruby laser for tattoo removal:

- Quick Warm-Up Time
- Large Spot Sizes
- Easy-to-Use Touchscreen Interface
- Optional MicroSpot Fractional Handpiece
- Superior Reliability

**Images Courtesy of Astanza Laser LLC.

ETERNITY^{TSR}

POWERFUL RUBY LASER FOR RESISTANT INKS

Technical Specifications

Wavelengths	694 nm
Indications	Tattoo Removal, Pigmented Lesions
Pulse Durations	40 ns
Fluence	Up to 20 J/cm ²
Repetition Rate	Up to 2 Hz
Spot Sizes	2.5, 4, 5, optional 6 mm
MicroSpot Fractional Handpiece	7.1mm x 7.1mm, 196 spots (14 x 14)
Electrical Requirement	110 V ± 10%, 50/60 Hz, 16 A



ASTANZA

SERVICE | TRAINING | MARKETING SUPPORT | TECHNOLOGY

Astanza is the innovator in tattoo removal technology, offering more cutting-edge tattoo removal solutions than any other aesthetic laser provider.

Focused on laser tattoo removal, Astanza knows what it takes to help practices succeed in this niche field. With the industry's fastest service support, in-depth training, and comprehensive marketing support, Astanza will be your partner every step of the way.

Corporate Headquarters

📍 1810 S Akard Street, Suite 500
Dallas, TX 75215

☎ (800) 364-9010

🌐 www.astanzalaser.com

